

HAPIShield

Cap To Help Prevent Hospital-Acquired Pressure Injuries

INSTRUCTIONS



SCAN THE QR CODE TO VISIT

Manufactured and Distributed by **INVENTIVE**
MEDICAL DEVICES

📍 3831 Turtle Creek Blvd., Unit 14F, Dallas Texas 75219. 📞 800-972-8855

HAPIShield - Description & Warnings

1. Description & Indications: Single-use patient care cap that attaches to O2 devices (nasal cannula, BiPAP masks, etc. that attach to the head or face) to assist with the prevention of HAPI (Hospital-Acquired Pressure Injuries).

2. Contraindications: HAPIShield Cap should only be used for securing hospital-grade oxygen supplies, i.e. nasal cannulas, HFNC, BiPAP masks, and other oxygen delivery devices with masks. HAPIShield strictly prohibits using the cap for any other purpose than listed above. Storage of the HAPIShield Cap should be between 50-80 degrees during storage and transport.

3. Warnings/Precautions:

- Single use only. Cleaning of devices is strictly prohibited by IMD or HAPIShield.
- If Velcro attachments become loose or break down, change the cap.
- Each patient may have different lengths of time to wear HAPIShield. Do not exceed 14 days of use for any single patient.
- External chemicals, gels, dyes, or stringent hair products may damage the cap.
- Any unauthorized use of HAPIShield will void the cap's warranty.
- Use only the recommended size of HAPIShield for the patient.
- HAPIShield will not be liable for any damages associated with oxygen delivery device malfunctions.
- HAPIShield will not be liable if the product is used past the indicated storage life shown on the product.
- Only expert clinical judgment should be used when sizing the HAPIShield Cap.

4. Adverse Reactions: The following adverse reactions may occur with HAPIShield:

- Skin abrasion from moving the product harshly against skin.
- Overtight or incorrect sizing of the cap.
- Patients who are not alert or awake must be checked regularly to assess sizing and skin integrity.
- Skin visualization is recommended every 2 to 4 hours to monitor for skin breakdown.

If any manufacturing defect occurs with HAPIShield during its use, please contact customer service. Regular clinical rounding must be performed with HAPIShield Cap to check skin condition and integrity.

HAPIShield - Instructions

Patient Set Up:

- 1) The cap will have the following:
 - A designated front and back of the cap.
 - Chin strap with Velcro for size adjustments.
 - Larger ear holes, with Velcro above the ears to secure tubing to the cap from O2 devices such as nasal cannulas, high-flow oxygen, oxygen masks, BiPAP, CPAP, and other oxygen delivery devices.
- 2) Place the cap on the head with "F" positioned at the middle of the forehead.
- 3) Lightly secure the chin strap.

For O2 Tubing - Ear Protection:

Place tubing on the patient as instructed by the device. Then, bring the tubing above the ear and secure it to the cap with Velcro.


For BiPAP/CPAP Protection:

Place the BiPAP or CPAP straps over the cap with a skin-protective barrier under the cap for high-risk patients.

Manufactured and Distributed by **INVENTIVE**
MEDICAL DEVICES

📍 3831 Turtle Creek Blvd., Unit 14F, Dallas Texas 75219. ☎ 800-972-8855

O2 TUBING

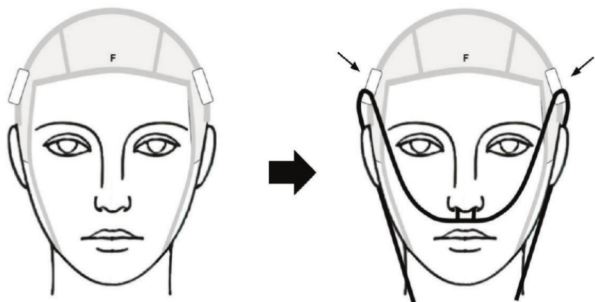

HIGH FLOW



BiPAP / CPAP

